



# PURELAB Option-S/R 7/15


The PURELAB Option provides pure water on demand with a wide range of flow-rates to suit your needs. Applications range from glassware washing and the feeding of ultra-pure water systems to cell culture and media preparation. The PURELAB Option-R provides additional bacteria and inorganic quality for your more critical applications.

- Quick and easy sanitization and replacement of consumables to reduce maintenance time
- Features 'Reverse Osmosis feed optimized' resin mixes to increase the ion exchange capacity of consumables and minimize running costs
- Delivers market leading microbial performance. In addition to simple sanitization and composite vent filtration, the PURELAB Option-R system is the only fully recirculating Type II system on the market, ensuring high specification water
- PURELAB Option systems are designed to be easy to access whether wall or bench mounted with a convenient dispense tap. They can be used with our wrap-around reservoirs to save space whilst optimizing purity



The only fully recirculating Type II system

### Process Flow PURELAB Option-R 7/15

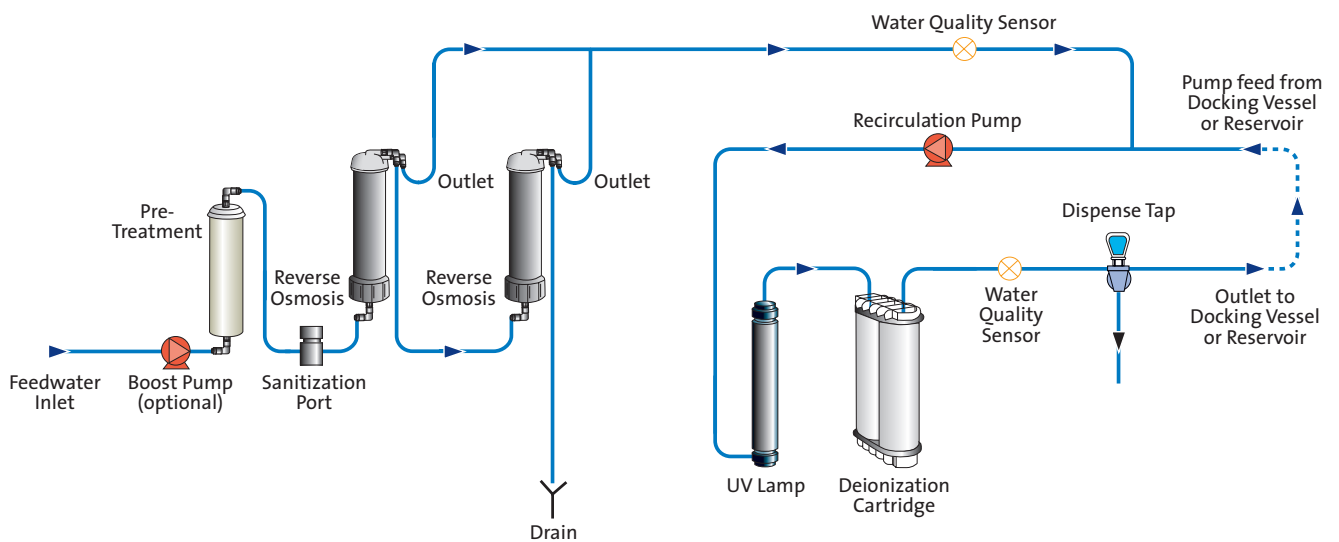

### Treated Water Specifications

Model	Option-S 7	Option-S 15	Option-R 7	Option-R 15
Make up rate <sup>1</sup>	7.5 l/hr	15 l/hr	7.5 l/hr	15 l/hr
Dispense rate from tap			1.0 l/min-nominal (less with POU filter)	1.0 l/min-nominal (less with POU filter)
Daily output (nominal max) <sup>1</sup>	180 l/24 hour day	360 l/24 hour day	180 l/24 hour day	360 l/24 hour day
Output reverse pressure (max) <sup>2</sup>	0.1 bar (1 psi)	0.1 bar (1 psi)	0.1 bar (1 psi)	0.1 bar (1 psi)
<b>Purity:</b>				
Inorganics@ 25°C	1 to >10 MΩ-cm	1 to >10 MΩ-cm	10 to >15 MΩ-cm	10 to >15 MΩ-cm
Total organic carbon (TOC)	<30 ppb	<30 ppb	<20 ppb	<20 ppb
Bacteria <sup>2</sup>	n/a	n/a	<1 CFU/ml	<1 CFU/ml
pH	Effectively neutral	Effectively neutral	Effectively neutral	Effectively neutral
Particles	n/a	n/a	Optional 0.2µm POU filter	Optional 0.2µm POU filter

<sup>1</sup> Standard conditions are 4 bar inlet pressure, 0 bar back pressure at 15 degrees centigrade, fed with potable water and a clean pre-treatment cartridge.

<sup>2</sup> Subject to suitable feedwater

### Dimensions and weights

Height	460mm (18.1in)	460mm (18.1in)	460mm (18.1in)	460mm (18.1in)
Width	410mm (16.1in)	410mm (16.1in)	550mm (21.7in)	550mm (21.7in)
Depth	270mm (10.6in)	270mm (10.6in)	270mm (10.6in)	270mm (10.6in)
Weight with internal boost pump	16kg (35lb)	16.5kg (36lb)	20kg (44lb)	21kg (46lb)
Weight without internal boost pump	13.5kg (30lb)	14.5kg (32lb)	18kg (40lb)	19kg (42lb)

### Feedwater Requirements

Source Quality	Potable mains water supply	Potable mains water supply	Potable mains water supply	Potable mains water supply
Fouling index - maximum	10	10	10	10
Total dissolved solids - maximum	1400µS/cm	1400µS/cm	1400µS/cm	1400µS/cm
Free chlorine - maximum	0.5 ppm	0.5 ppm	0.5 ppm	0.5 ppm
Heavy metals - maximum	0.05 ppm	0.05 ppm	0.05 ppm	0.05 ppm
Silica - maximum	30 ppm	30 ppm	30 ppm	30 ppm
Temperature	1 - 35°C	1 - 35°C	1 - 35°C	1 - 35°C
Flowrate (maximum requirement)	78 l/hr	85 l/hr	78 l/hr	85 l/hr
Drain requirements (gravity fall with air gap). Maximum during service	70 l/hr			
<b>Feedwater Pressure</b>				
Maximum - without internal boost pump	6.0 bar (90 psi)	6.0 bar (90 psi)	6.0 bar (90 psi)	6.0 bar (90 psi)
Minimum - without internal boost pump	4.0 bar (60 psi)	4.0 bar (60 psi)	4.0 bar (60 psi)	4.0 bar (60 psi)
Maximum - with internal boost pump	2.0 bar (30 psi)	2.0 bar (30 psi)	2.0 bar (30 psi)	2.0 bar (30 psi)
Minimum - with internal boost pump	Flooded Suction	Flooded Suction	Flooded Suction	Flooded Suction

### Electrical Requirements

Mains input	100-240V ac, 50-60Hz	100-240V ac, 50-60Hz	100-240V ac, 50-60Hz	100-240V ac, 50-60Hz
System voltage	24V dc	24V dc	24V dc	24V dc
Power consumption with boost pump	43VA	43VA	80VA	80VA
Power consumption without boost pump	19VA	19VA	50VA	50VA
Fuses	2 x T3.15 Amp	2 x T3.15 Amp	2 x T6.3 Amp	2 x T6.3 Amp
Reservoir level connection	Jack Plug 3.5mm	Jack Plug 3.5mm	Jack Plug 3.5mm	Jack Plug 3.5mm
Noise level	<45dBA	<45dBA	<45dBA	<45dBA

### ELGA LabWater

Tel: +44 (0) 1494 887500 Fax: +44 (0) 1494 887505

Email: [info@elgalabwater.com](mailto:info@elgalabwater.com) Website: [www.elgalabwater.com](http://www.elgalabwater.com)

ELGA® is the global laboratory water brand name of Veolia Water. VWS (UK) Ltd. Registered in England & Wales No. 327847 ©Copyright 2008 ELGA LabWater/VWS (UK) Ltd. All rights reserved. As part of our policy of continual improvement we reserve the right to alter the specifications given in this datasheet.

LITR 38753-01